

02.


EXTRACORONAL ATTACHMENTS

CEKA 334
ENGLISH



CEKA
ATTACHMENTS
PRECI-LINE


**PLENTY OF SPACE
FOR AN ATTACHMENT?**

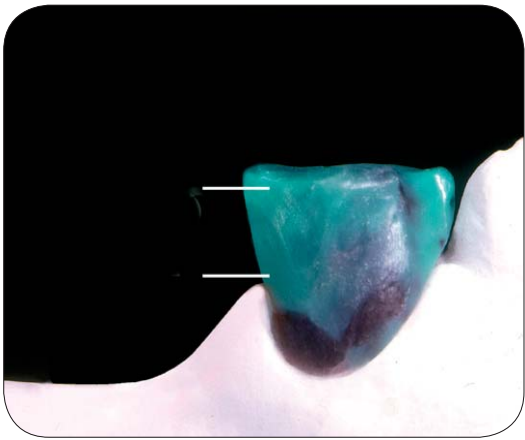

**OF COURSE,
WITH CEKA REVAX!**




FAMILIAR COMPONENTS

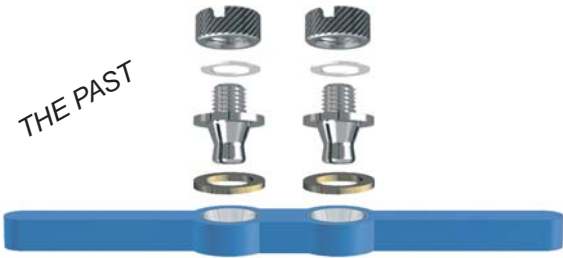

OL 334 NX

2 attachments/package

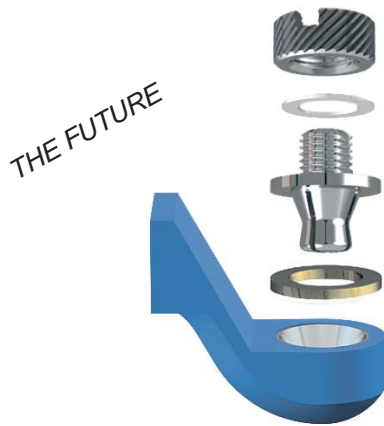

OL 0234 NPR

1 attachment/package

TRY IT!



Press the cast-on precision female with the colour-coded side down into the keeper.



Reduce the plastic keeper at the buccal. This creates a natural transition to the crown and more space for aesthetics.



Position the profile on the crest with passive contact. Cast in a non-precious alloy.



The attachment female provides adequate space for good contour and aesthetics for a natural appearing removable prosthesis.



CEKA AXIAL
PRECI-CLIX AXIAL
PRECI-CLIX RADICULAR
PRECI-BALL

CEKA EXTRACORONAL
CEKA ID
PRECI-VERTIX
PRECI-VERTIX AT
PRECI-CLIX EC
PRECI-SAGIX
PRECI-52
PRECI-TUBIX/MORTIX

PRECI-BAR
PRECI-HORIX
PRECI-CLIP

PRECI-PROFILE
PRECI-POST
CEKA SOL
CEKALLOY
CEKA SITE
CEKA BOND
3C-BOND
PERMA-RET
PRECI-SEP
PLASTICWAX
EXPANDO

01.
AXIAL
ATTACHMENTS

02.
EXTRACORONAL
ATTACHMENTS

03.
PROFILE
ATTACHMENTS

SPECIALTIES

WWW.CKPL.EU

CEKA
ATTACHMENTS
PRECI-LINE



Prospectus n° 263.E-ed.B

© Copyright ALPHADENT NV 10/2011 - BTW BE 0414.789.321 RPR Kortrijk

ALPHADENT NV • Textielstraat 24 8790 Waregem, Belgium

T+32 56 629 100 - **F**+32 56 629 101 - **E**info@ckpl.eu

